Implantlink Instructions

implantlink® semi Forte

GB
Instructions for use

Resin-based, temporary luting cement with increased adhesion for implant-supported restorations

Dual-curing, semi-permanent implant cement with increased adhesion and compressive strength. Specially formulated for customised abutments that have a reduced adhesive area or are particularly small and low. Extremely low film thickness thanks to very low displacement resistance. The highly polymerised, ductile resin structure of the material gives it high compressive strength with a high marginal seal. Residues can be easily removed in large pieces in the plastic phase without leaving undesired particle residues. The attached restorations can be removed without damage thanks to the perfectly balanced adhesion properties of the material. For all combinations of material. Eugenol-free, antibacterial and neutral in odour and taste.

Fig. 1

1. Preparation and cleaning of implant abutments and inner surfaces of restorations

Carefully check fit, occlusion and proximal contact points. Complete all finishing procedures. Degrease, clean (e.g. alcoholic solutions) and dry prepared abutments. Make sure the interior surfaces of the superstructure are clean and dry. (Fig. 1).

Fig. 2

2. Mixing and applying

Bleed the syringe before installation of the automix cannula until material emerges uniformly from both openings (Fig. 2). Insert the plunger in the body of the cartridge. Remove the cap from the cartridge by turning clockwise and attach mixing cannula. Pay attention to the guides on cannula and cartridge. Lock in place by turning clockwise (Fig. 3). Again, bleed small amount to ensure uniform and even material flow. Due to finely tuned dispensing pressure, dosing can be performed individually and precisely. Dispense the material with uniform pressure. After use, leave mixing cannula on cartridge until next application to keep cement from drying out.

Fig. 3

3. Inserting the restoration

Apply a thin layer of **implantlink® semi Forte** directly into the superstructure (Fig. 4). Carefully position the restoration on the implant abutment with slight pressure (Fig. 5); then press down firmly. Remove excess material during the gel phase, after around 2 – 3 minutes, with a suitable instrument (Fig. 6).

Fig. 4

4. Working time / light curing

The working time at room temperature is around 80 seconds. The cement is cured at body temperature in about 5 – 6 minutes. Direct exposure to light accelerates setting time. Light curing (marginal region, excess material) is achieved with approx. 20 seconds per surface (Fig. 7). When **implantlink® semi Forte** is used with translucent materials (e.g. thin layers of ceramic provisional crowns and bridges), increase exposure to light through the material to 30 - 40 seconds per surface.

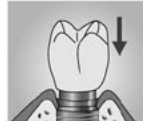
Fig. 5

Precautions

- Do not rough or sand-blast abutment surface.
- Do not damage the underlying titanium surface when removing excess material. If necessary, use plastic instruments (implant scalers).
- Leave no cement residues in the mouth or gingival sulcus.
- Use no adhesion-reducing additives. They could uncontrollably render the retentive qualities of **implantlink® semi Forte**.
- Do not manually remove or mix cement components.
- Avoid contact of unmixed components with skin or mucosa.
- In very rare cases, **implantlink® semi Forte** may cause sensitization. Should such reaction occur, discontinue the use of this material.
- Observe the safety data sheet.
- Do not use **implantlink® semi Forte** for the luting of hybrid restorations with natural teeth and implant abutments.
- Only for intended use by qualified personnel.

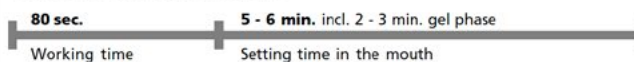
Fig. 6

Fig. 7

Autopolymerisation:

Light-curing:

DETAX

GmbH & Co. KG
Carl-Zeiss-Str. 4 · 76275 Ettlingen/Germany
www.detax.de · post@detax.de

more information
www.detax.de/en
implantlink/clip

Made in Germany

02/2011

Indication for Use:

- Temporary luting of implant-supported temporary and permanent restorations on customised, reduced and small implant abutments

Technical data:

- Mixed volume:** 5 ml (mini-mix)
- Mixing ratio:** 4:1
- Colour code:** Base: white-opaque Catalyst: semi-transparent
- Mixing time:** None: automix (mini-mix)
- Working time:** approx. 80 sec.
- Gel phase:** 2 - 3 min. * (for removal of excess material during placement)
- Setting time in the mouth:** 5 - 6 min. *
- Setting time with light-curing:** approx. 20 sec. per surface; 30-40 sec. through translucent materials
- Film thickness:** < 10 µm
- Compression strength:** > 100 MPa
- Application temperature:** At 23 °C ± 2 °C / 73 °F ± 4 °F; 50 ± 5 % rel. humidity
- Storage:** below 25° C / 77° F

* from beginning of mixing at 23 °C ± 2 °C / 73 °F ± 4 °F; 50 ± 5 % rel. humidity. Higher temperatures accelerate, lower temperatures retard a. m. times.

Ordering information:

implantlink® semi Forte
Standard packing **03371**
cartridge of 5 ml mini-mix 4:1
10 mixing cannulas, brown

implantlink® semi Classic
Standard packing **03092**
cartridge of 5 ml mini-mix 4:1
10 mixing cannulas, brown

Mixing cannulas 4:1 **02591**
brown, 25 pcs.

monoprint supra
Transfer silicone
Standard packing **02525**
2 x 50 ml cartridges
6 mixing cannulas, green

Jumbo pack 5:1 03377
maxi cartridge of 300 ml
base + 60 ml catalyst
6 dynamic mixers, yellow
1 rotating ring, yellow

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